




Home care masks and accessories

A solution for every patient

PHILIPS


RESPIRONICS

Which mask for which patient?

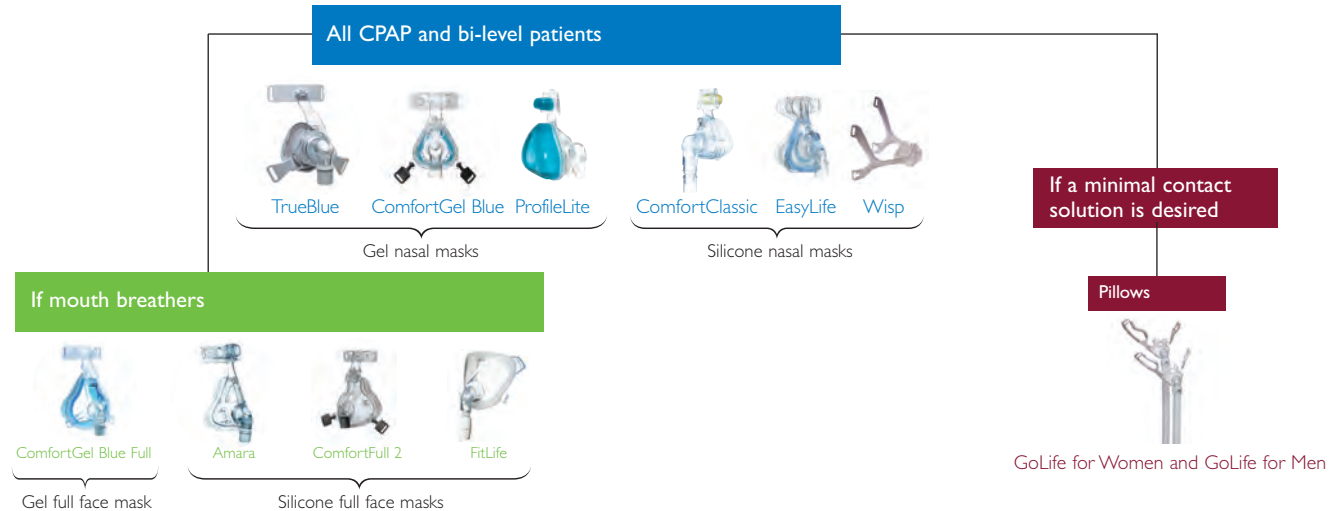

Nasal masks

Full face masks

Minimal contact masks

Accessories

Mask fitting

Devices

Patient interface line overview

TrueBlue



ComfortGel Blue



ProfileLite



ComfortClassic EasyLife



Wisp



Nasal masks

Sizes	P, S, M, MW, L	P, S, M, L FitPack (S/M) DuoPack (P,S,M,L)	P, S, SM, M, ML, L, LN, Youth, Child	S, M	P, S, M, MW, L FitPacks (S/M, MW/L) DuoPacks (S, M, MW, L)	P, S/M, L
SE version (no exhalation port)	Yes	Yes	Yes	No	Yes	Yes
Replaceable cushion	Yes	Yes	No	No	Yes	Yes

Amara



ComfortGel Blue Full



ComfortFull 2



FitLife



GoLife



Full face masks

Sizes	P, S, M, L	S, M, L, XL	S, M, L	S, L, XL
SE version (no exhalation port)	Yes	Yes	Yes	Yes
Replaceable cushion	Yes	Yes	No	No

Minimal contact mask

Female (P, S, M) Male (S, M, L)
No
Yes

Nasal masks

Full face masks

Minimal contact masks

Accessories

Mask fitting


Devices

TrueBlue


The next generation gel nasal mask with AutoSeal technology



Without exhalation port

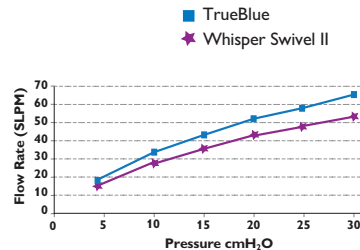

Size	P	S	M	MW	L
With exhalation	1071838	1071839	1071840	1071841	1071842
Without exhalation	1071843	1071844	1071845	1071847	1071848

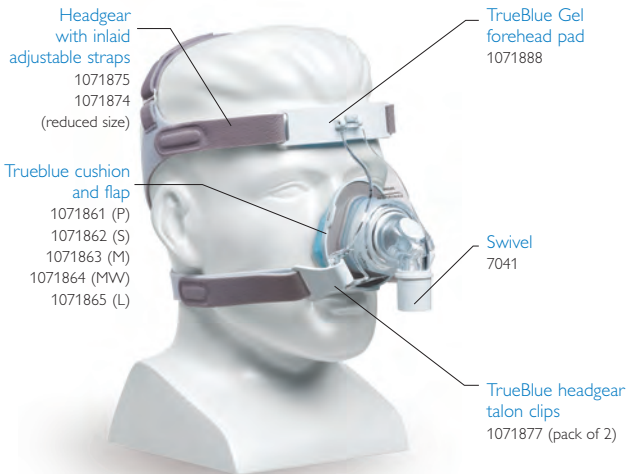


Single use sizing gauge: 1071878

ComfortGel Blue

Our gel technology for an optimized comfort and performance



Without exhalation port

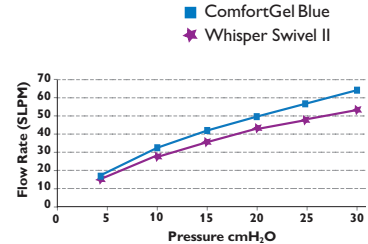

Size	P	S	M	L
FitPack S/M	1070076			
DuoPack	1070075	1070074	1070073	1070072
With exhalation	1070067	1070066	1070065	1070064
Without exhalation	1070071	1070070	1070069	1070068




ProfileLite

The original gel mask without silicone and fewer adjustment features




Without exhalation port

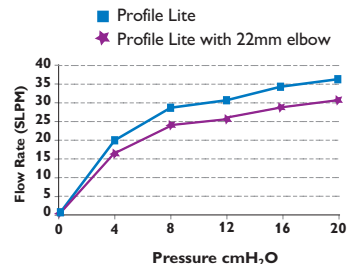


Size	P	S	MS	M	MW	L	LN	Youth	Child
With exhalation	1004113	1004114	1004115	1004116	1004117	1004118	1004119	1054707	1008829
Without exhalation	1006314	1006315	1006316	1006317	1006318	1006319	1006320		

ComfortClassic

Our simple silicone mask with very few adjustments

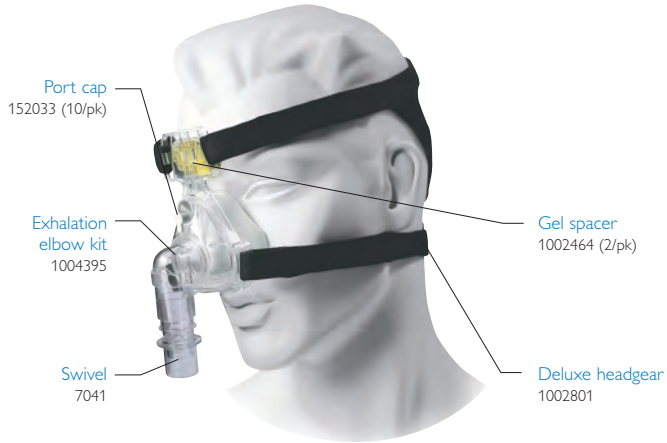

Comfort Classic sizing gauge: 1012041

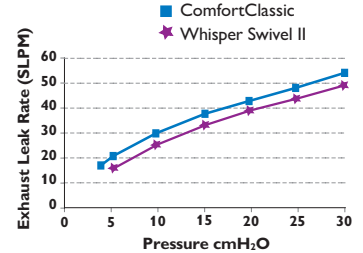

Size	S	M
With exhalation	1007967	1007968


Nasal masks

EasyLife


Our self-adjusting silicone nasal mask with AutoSeal technology



Without exhalation port

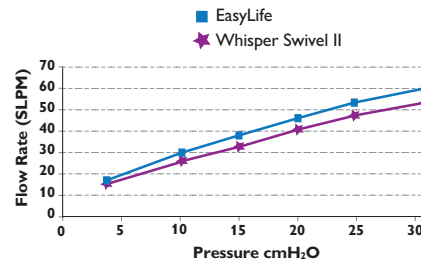



Size	P	S	M	MW	L
With exhalation	1050020	1050021	1050022	1050023	1050024
Without exhalation	1050030	1050031	1050032	1050033	1050034
DuoPack (extra cushion)	1050025	1050026	1050027	1050028	1050029
FitPack (mask with headgear and cushions S and M)	1050088				
FitPack (mask with headgear and cushions MW and L)	1050089				

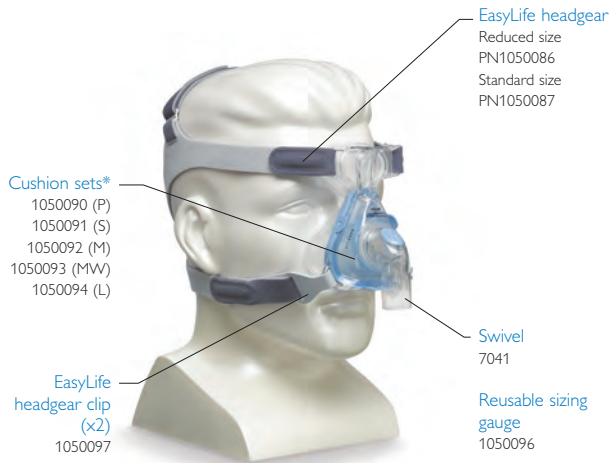

*Cushion sets include both inner and outer cushions of the same size

Wisp

Performance has a new style

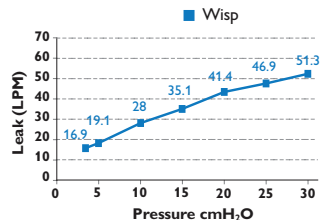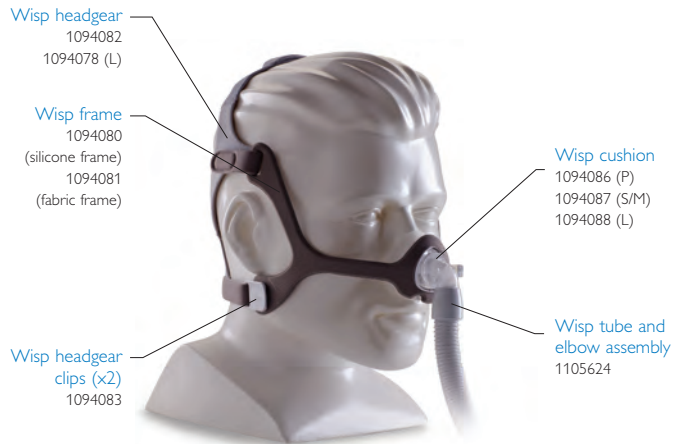

With exhalation port	Pack with cushion set*
Wisp mask with silicone frame	1094060
Wisp mask with fabric frame	1094061
Without exhalation port	
Wisp mask with silicone frame and without exhalation	1094062



*Each package contains all three cushion sizes, which fit on the same frame.




Nasal masks

Amara

Performance, appeal and simplicity at its best.

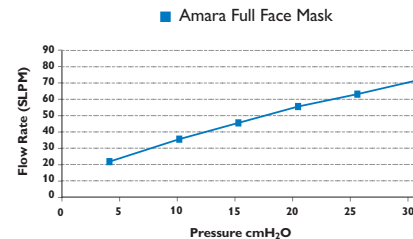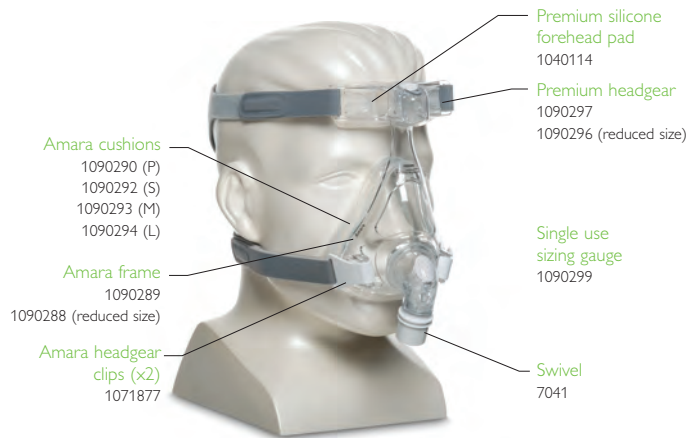

Size	P-RS	S-RS	S	M-RS	M	L
With exhalation	1090220	1090221	1090222	1090225	1090223	1090224
Without exhalation	1090230	1090231	1090232	1090235	1090233	1090234

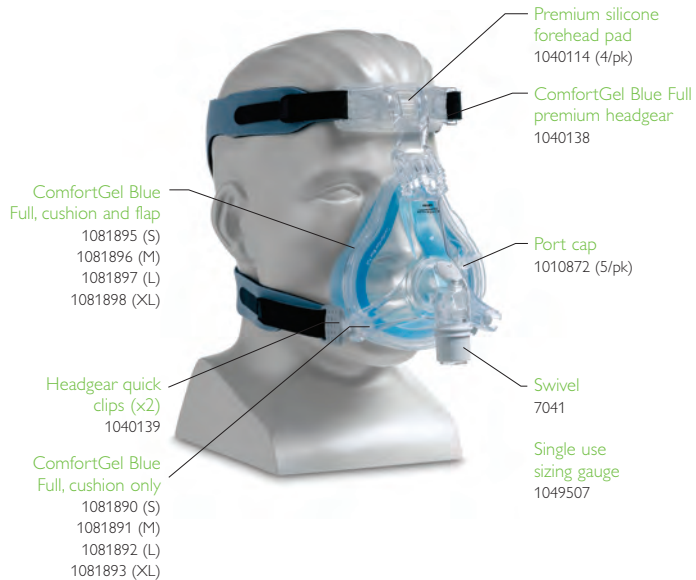


Full face masks

ComfortGel Blue Full


Our gel technology now available for mouth breathers.



Without exhalation port

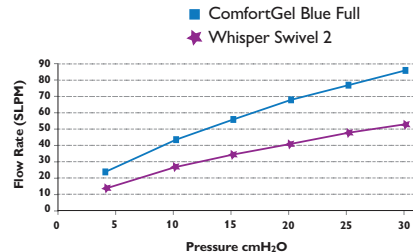

Size	S	M	L	XL
With exhalation	1081820	1081821	1081822	1081823
Without exhalation	1081830	1081831	1081832	1081833

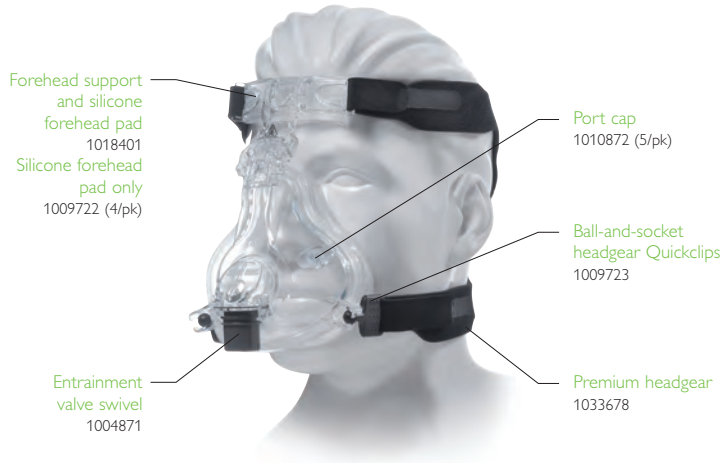


ComfortFull 2

The basic silicone full face mask for mouth breathers



Without exhalation port

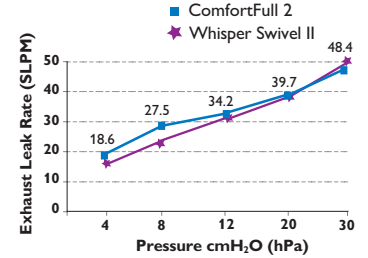


Size	S	M	L
With exhalation	1004881	1004873	1004951
Without exhalation	1023747	1024261	1024262

FitLife

A total face mask that provides comfort for hard-to-fit patients



Without exhalation port

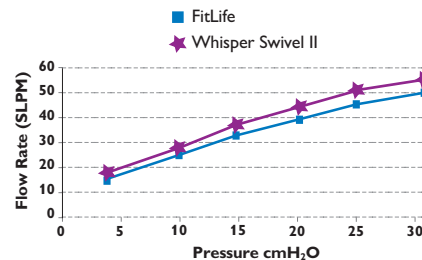


Size	S	L	XL
With exhalation	1060803	1060804	1089994
Without exhalation	1061711	1061712	1089995

Full face masks



GoLife

The unique fit each one deserves

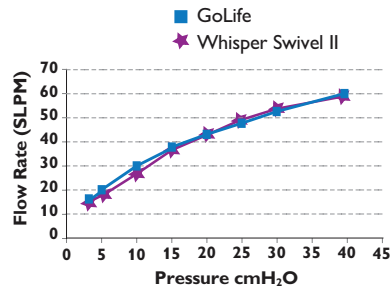

Size	Pack with cushion set
For Women (P, S, M)	1082580
1 frame, 1 headgear, 1 set of cushions	
For Men (S, M, L)	1073331
1 frame (+ 1 additional smaller frame), 1 headgear, 1 set of cushions	

Minimal contact masks

Accessories

Headgear

Part no.	Headgears
1033678	Premium headgear with EZ Peel tabs (all headgear fit any nasal or fullface mask – except for FullLife and FitLife)
1059257	Premium smaller size headgear
1040138	Blue Premium headgear with EZ Peel tabs (alternative headgear for ComfortGel Full mask)
1002801	Deluxe headgear (alternative headgear for ProfileLite and ComfortClassic masks)
1018755	Deluxe smaller size headgear
1002285	SimpleStrap headgear
1050087	EasyLife headgear
1050086	EasyLife headgear reduced size
1071875	TrueBlue headgear
1071874	TrueBlue headgear reduced size
1060807	FitLife headgear, S
1060808	FitLife headgear, L
1089990	FitLife headgear, XL
1094082	Wisp headgear
1094078	Wisp headgear, L



1033678



1040138



1002801



1050095



1060808

Softcaps

Part no.	Softcaps (only fit to ProfileLite, ComfortClassic, and Contour masks)	
Mesh softcap		
1001632	Mesh softcap - Child	P
1002806	Mesh softcap - Medium size	M
1002807	Mesh softcap - Large size	L
1001634	Mesh softcap - X-Large size	XL
White Lycra softcap		
1001624	White softcap - Child	P
1001625	White softcap - Medium size	M
1001626	White softcap - Large size	L
1001627	White softcap - X-Large size	XL
Blue Lycra softcap		
1001628	Blue softcap - Child	P
1001629	Blue softcap - Medium size	M
1001630	Blue softcap - Large size	L
1001631	Blue softcap - X-Large size	XL



Mesh softcap



White Lycra softcap



Blue Lycra softcap

Chin straps

Part number	Chin straps
1012911	Premium chin strap
1006085	Deluxe chin strap
302175	Chin strap



1012911



1006085



302175

Reusable circuits

Part number	Reusable circuits
622038	Reusable flexible tubing, grey (1,83 m)
1032907	Performance tubing 22 (22mm)
1074400	Performance tubing 15 (15mm)
1049270	Tubing insulator (for rainout and noise reduction)
1007800	High-temperature tubing, beige (1,83 m)
1006121	Reusable patient tubing with pressure line (1,83 m)
7212	Reusable grey tubing (Long: 3 m)



622038



1032907



1049270



1007800



1006121

Reusable valves and swivels

Part number	Reusable valves and swivels
302312	Plateau exhalation valve (PEV)
332113	Whisper Swivel 2
622148	Outlet swivel with pressure tubing
302418	Pressure valve



302312



332113



622148



302418

Other accessories

Part number	Other accessories
312710	O ₂ enrichment port
1014047	Bacteria filter
FC06566	Mouthpiece, angled, 22 mm
1004524	Mouthpiece, angled, 15 mm
452040	NG tube sealing pads (10/Pk)




312710



1014047



FC06566

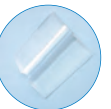

452040


3 steps to successful mask fitting

Size selection – step 1

Nasal mask sizing

Landmarks

- a. Sides of nostrils
- b. Bridge of nose zone
- c. Below the nose tip, above the lip

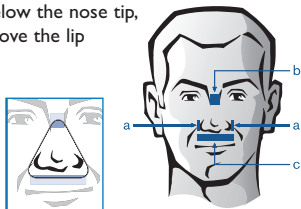

← Use the sizing gauge to select the right size →



or



Applicable to TrueBlue only




Full face mask sizing

Landmarks

- a. Sides of mouth
- b. Bridge of nose zone
- c. Below the lower lip

Make sure the mouth is slightly open when checking size and the top of the cushion is not too high on the bridge of the nose area.

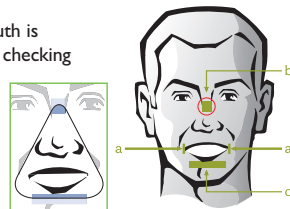


For individual mask fitting steps, please refer to the user instruction guide of your mask.

Mask fitting

Mask fitting – step 2

Common mistake
StabilitySelector adjusted
too high / bottom strap
too tight = leaks into
eyes, top lip pressure /
irritation

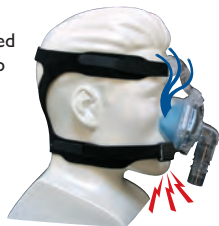


Common mistake
StabilitySelector adjusted
too low / bottom strap
too tight = leaks around
the mouth, forehead and
bridge nose irritation /
pressure points




How to adjust?

1. Click and slide the adjuster downwards.
2. Loosen the bottom strap of headgear.

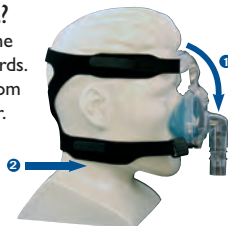


How to adjust?

1. Click and slide the adjuster upwards.
2. Loosen top strap of headgear.

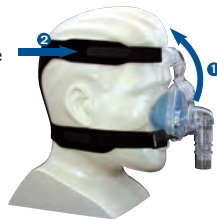

One click on the StabilitySelector may make the difference.
Always move the StabilitySelector in the opposite direction
of the leak.

OK

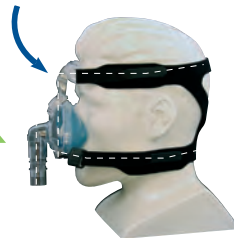

Correct position = mask parallel to the face
Do not overtighten the headgear.

For individual mask fitting steps, please refer to the user instruction guide of your mask.

Final adjustments and cleaning instructions – step 3

Check the adjustments of the mask in lying position with the right pressure levels.



Cleaning instructions

Mask and tubing (wash daily)

- Hand wash the mask in warm water with mild dishwashing detergent.
- Rinse thoroughly and air dry completely before use.
- Wash face carefully before each use.
- Inspect the mask after cleaning; replace the mask and nasal interface cushion if damaged or worn.
- Rinse the tubing under shower and hang to dry.

Headgear (wash weekly)

- Handwash the headgear with mild detergent and warm water.
- Rinse thoroughly. Pat with a towel to absorb excess moisture, then hang to dry.

⚠ DO NOT use bleach, alcohol, any strong household cleaners, or cleaning solutions with alcohol, conditioners or moisturizers.

Additional fitting tips

- Check that the headgear is not overtightened: allow one finger/hand space between the face and each of the straps.
- Reposition the mask, if necessary, by gently pulling the cushion away from the nose and then back to the initial fit.





















For individual mask fitting steps, please refer to the user instruction guide of your mask.


Mask fitting

Which device for which patient?

device for which patient?

Sleep therapy				Home ventilation		
   				  		
          						
System One REMstar Pro				BiPAP S/T		
System One REMstar Auto				BiPAP A30 and A40		
System One BiPAP Auto				Trilogy100		
System One BiPAP autoSV Advanced						
Target patients	OSA			Complex breathing patterns such as: CSR, CSA, mixed apnea, and complex apneas		
	OSA			Non-dependent chronic respiratory insufficiency patients (adult and children > 7 yrs old and 18 kg)		
Modes	Fixed CPAP (CPAP check)			Non-dependent chronic respiratory insufficiency patients (adult and children > 10 kg)		
	Fixed CPAP Auto CPAP (CPAP Check and auto-trial)			Dependent chronic respiratory insufficiency patients (adult and children > 5 kg)		
	Fixed BiPAP Auto BiPAP			Setup possible as: Auto CPAP, BiPAP S and S/T, Auto CPAP + Servo Ventilation, BiPAP + Servo Ventilation		
				• BiPAP S/T: CPAP, S, S/T		
				• BiPAP A30: CPAP, S, S/T, PC, T • BiPAP A40 (NIV and IV capabilities, integrated battery*): CPAP, S, S/T, PC, T, AVAPS-AE		
			• Invasive and non-invasive capabilities • Volume modes: AC, CV, SIMV • Pressure modes: CPAP, S, S/T, PC, T, PC-SIMV			

 Humidifier available in option

 Heated hose available in option

*Battery module may be optional in some markets.

Devices

Notes

Not all products are available in all countries.
Please check with your local Philips Respironics
representative for availability.

Respironics and the name of Respironics products contained in this document are trademarks of Respironics, Inc. and its affiliates. All rights reserved.



©2013 Koninklijke Philips Electronics N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.

Broudy AE 5/8/13 MCI 4105504 PN 1107381

Philips Healthcare is part of Royal Philips Electronics

Philips Respironics International Headquarters
+33 1 47 28 30 82
Philips Respironics Asia Pacific
+65 6882 5282
Philips Respironics Australia
+61 (2) 9947 0440
Philips Respironics United Kingdom
+44 800 1300 845

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service
+1 724 387 4000
800 345 6443 (toll free, US only)
www.philips.com/respironics